


抗てんかん薬 Antiepileptics

Key words : 抗てんかん薬 医薬品
Column : USP category: L1


[Analytical conditions]

Column : L-column ODS (C18, 5 μm, 12 nm), 4.6 mm I.D. × 250 mm L.; Cat. No. 622080
Eluent : CH₃CN/20 mM Phosphate buffer pH 6.8 (32/68)
Flow rate : 1 mL/min
Temperature : Room temperature
Detection : UV 220 nm
Injection volume : 25 μL
System : -

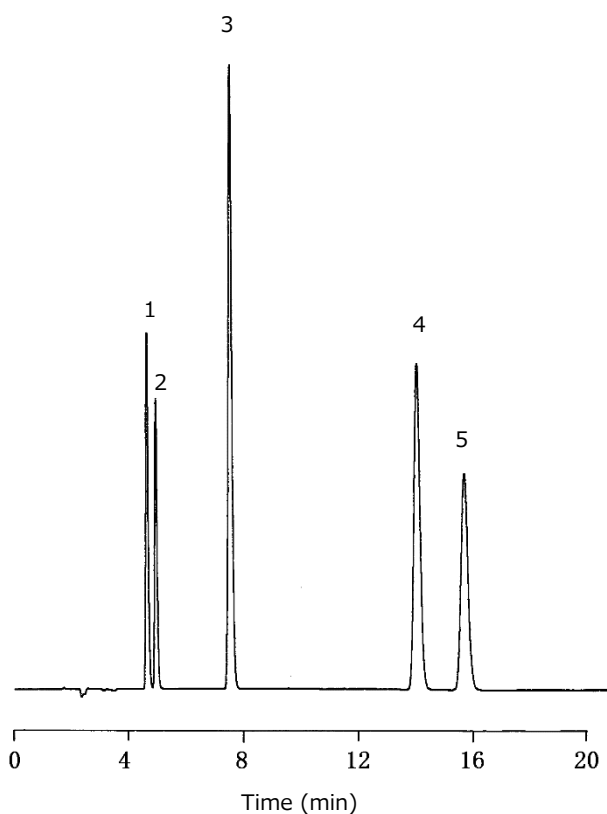


Sample:


1. Primidone (10 mg/L)
プリミドン

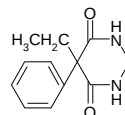


2. Ethosuximide (100 mg/L)
エトスクシミド

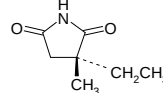

3. Phenobarbital sodium salt (20 mg/L)
フェノバルビタールナトリウム

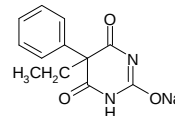

4. Phenytoin (10 mg/L)
フェニトイン

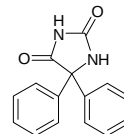

5. Carbamazepine (5 mg/L)
カルバマゼピン

